



# Therapy you can trust

**InnoSpire Deluxe** – compressor nebulizer system

InnoSpire Deluxe is a reliable and easy to use compressor nebulizer system suitable for nebulizing most aerosol medications prescribed for common respiratory diseases.

InnoSpire Deluxe is part of our InnoSpire family of compressor nebulizer systems that provide reliable drug delivery to patients.

Suitable for professional or home use, InnoSpire Deluxe has been designed to easily manage prolonged use in a clinical environment as well as frequent use in the home.

#### **Proven nebulizer technology**

Coupled either with the SideStream Reusable or SideStream Plus handset, the InnoSpire Deluxe provides highly efficient aerosol therapy in fast treatment times.

#### **Key features**


- Matched with SideStream Reusable or SideStream Plus
- 7LPM driving flow
- Continuously rated for constant use
- Large integral storage
- Handle for easy transport

**PHILIPS**  
RESPIRONICS

# Reliable aerosol delivery

**SideStream Reusable and SideStream Plus have been designed to draw in more air, via an active venturi system, resulting in fast drug delivery and short treatment times.**

They also have a unique diamond jet design that helps to minimize wear from frequent use and aids consistent delivery of a high quality aerosol.

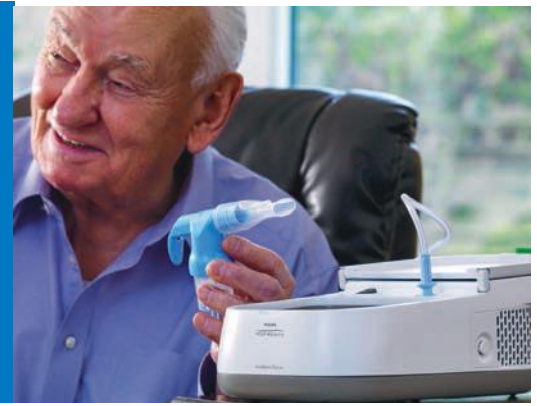

## SideStream Reusable

By utilizing a venturi, SideStream produces a flow rate of 16 l/min



## SideStream Plus

Equipped with a low resistance valve to deliver fast treatments with less waste



Specifications	Product	
<b>Performance data</b>	<b>SideStream Plus</b>	<b>SideStream Reusable</b>
MMAD	3.42µm	2.93µm
Fine particle fraction	67%	76.3%
Treatment time (2.5mL charge)	3.68 min	4.12 min
<b>Product specifications</b>	<b>InnoSpire Deluxe</b>	
Electricity supply	230V/50Hz, 1Amp	
Compressor max. pressure	360kPa (3.6bar)	
Average flow rate	7LPM @ 10 psig	
Max flow rate	9.8LPM	
Thermal overload protector	Thermal fuse; functioning temperature 150°C	
Size	38 x 17 x 11cm	
Weight	2.2kg	
Sound level	58 ± 3dBA	
Operating temperature/humidity	10°C to 40°C / 10% to 95% non-condensing	

Ordering information	Part no.
<b>Compressor nebulizer systems - including SideStream Plus kit</b>	
InnoSpire Deluxe – 230V/50Hz – EU plug	1110057
InnoSpire Deluxe – 230V/50Hz – UK plug	1110064
<b>Compressor nebulizer systems - including SideStream Reusable kit</b>	
InnoSpire Deluxe – 230V/50Hz – EU plug	1112278
InnoSpire Deluxe – 230V/50Hz – UK plug	1112279
InnoSpire Deluxe – 220V/60Hz – UK plug	1110061
InnoSpire Deluxe – 230V/50Hz – EU plug	1110058
InnoSpire Deluxe – 230V/50Hz – EU plug	1110065
<b>Accessories/replacement parts</b>	
SideStream Plus reusable kit	1092001
SideStream reusable year pack	2010A
SideStream reusable kit with adult mask	1223A
SideStream reusable kit with pediatric mask	1224A
SideStream reusable kit with mouthpiece	1225A
Air inlet filters (pack of 4)	1113525

Please visit [www.philips.com/innospire](http://www.philips.com/innospire) for more information



© 2013 Koninklijke Philips N.V.  
All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

Broudy NB 12/17/13 MCI 4105934 PN 1113989

Philips Healthcare is part of Royal Philips

[www.philips.com/healthcare](http://www.philips.com/healthcare)  
[healthcare@philips.com](mailto:healthcare@philips.com)